- The installation location must be within the reception range of the Carephone. There must be good radio signal reception. See Initial Operation.
- The detector is only suitable for living spaces.
- The detector must not be directed toward window surfaces.
- Do not install in direct sunlight.
- Do not direct toward objects that experience significant changes in temperature.
- Do not direct toward rotating machines.
- Installation on metal surfaces can impair the range of the radio signal.

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**!** Installation and initial operation should only be carried out by trained service personnel.

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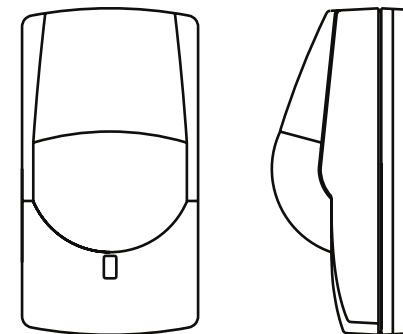
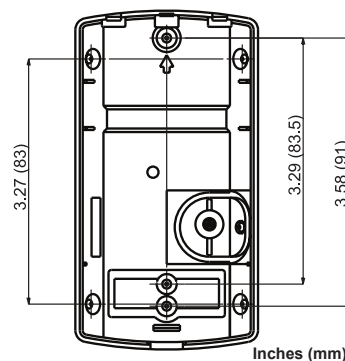
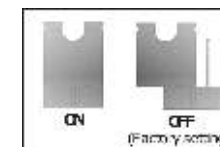
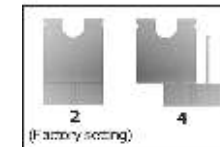
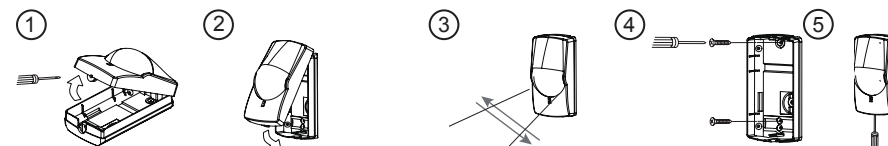
**EN Operating Instructions**

The Wireless Motion Detector (passive infrared detector) is used to monitor a room and to wirelessly transmit presence information to the Carephone. Unlike triggering via the daily button, there is no audio announcement if triggering occurs via the detector.

**Jumper block functions**

**LED ON/OFF** The detector LED indicator can be switched either "ON" or "OFF"

**PULSE COUNT 2/4** The detection mode can be selected depending on the environmental conditions of installation. Select "2" for normal applications or select "4" for use in areas where there may be movement of small animals (<40cm).

**WALK TEST** is used to control the detection coverage. The LED will illuminate at each detection.

**Initial Operation**

Put the Carephone in learning mode to record the wireless motion detector WNX-40PI-T.

- (1) Open the detector housing using a screwdriver.
- (2) Place the jumper block indicated "WALK TEST" (closed). Then, install the battery in the battery holder.
- (3) Pass your hand in front of the detector; the LED will illuminate for 3 seconds and the event is transmitted to the Carephone.

**This mode is only possible for 2 minutes after the first battery insertion. After this time, as soon as a movement has been reported to the Carephone, no further events will be transmitted within the next 12 minutes.**

**Installation**

- (4) Fix the bottom of the housing on a wall according to the information on the first page.

Remove the jumper block "WALK TEST" (normal position = battery saving mode) and close the housing detector.

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