



VeriPass



VeriPass
Klik hier als u tekst wilt invoeren.

11/09/2013

Copyright

This document may not be copied in whole or in part or otherwise reproduced without prior written consent from IndigoCare Europe or IndigoCare USA, except where specifically permitted under international copyright law.

Disclaimer

The information in this document is subject to change without notice. IndigoCare Europe or IndigoCare USA assumes no responsibility for inaccuracies or omissions and specifically disclaims any liabilities, losses, or risks, personal or otherwise, incurred as a consequence, directly or indirectly, of the use or application of any of the contents of this document. For the latest documentation, contact IndigoCare Europe or Indigocare USA.

This publication may contain examples of screen captures and reports or database used in daily operations. Examples may include fictitious names of individuals and companies. Any similarity to names and addresses of actual businesses or persons is entirely coincidental.

All examples of patient information in this document are fictitious. Any resemblance to a real patient or facility is purely coincidental. The owners and users of this product are solely responsible for complying with all applicable patient information laws. The users, by their use of this product, agree to indemnify the manufacturer and/or seller of this product against all claims, litigation, and suits filed for patient information violations.

Intended use

Use this product only for the purpose it was designed for; refer to the data sheet and user documentation for details. For the latest product information, contact Indigocare Europe or IndigoCare USA.

This guide is a working and supporting tool for certified engineers for installing and maintaining iCall installations.

Contact information

IndigoCare USA
1001 South Dairy Ashford, Suite 100
Houston, TX 77077
USA

T: +1 713.595.6698
F: +1 713.595.6699
Email: usa@indigocare.com
www.indigocare.com

IndigoCare Europe
Schoebroekstraat 48
3583 Paal
Belgium

T: +32 11 247090
F: +32 11 247099
Email: Europe@indigocare.com
www.indigocare.com


Content

- 1. General purpose5
- 2. Login5
- 3. Menu6
- 3.1. Scan the network for access control units.....7
- 3.2. Manage groups8
- 3.3. Manage users..... 11
- 3.4. Manage rooms..... 13
- 3.5. Manage access rights 15
- 3.6. Upload all settings to the system 16
- 3.7. Settings screen 17

1. General purpose

This document describes the implementation, programming and testing of the Access Control software. The goal is to send access control rights to the appropriate iCall IP units, to handle the access control in a correct manner.

2. Login

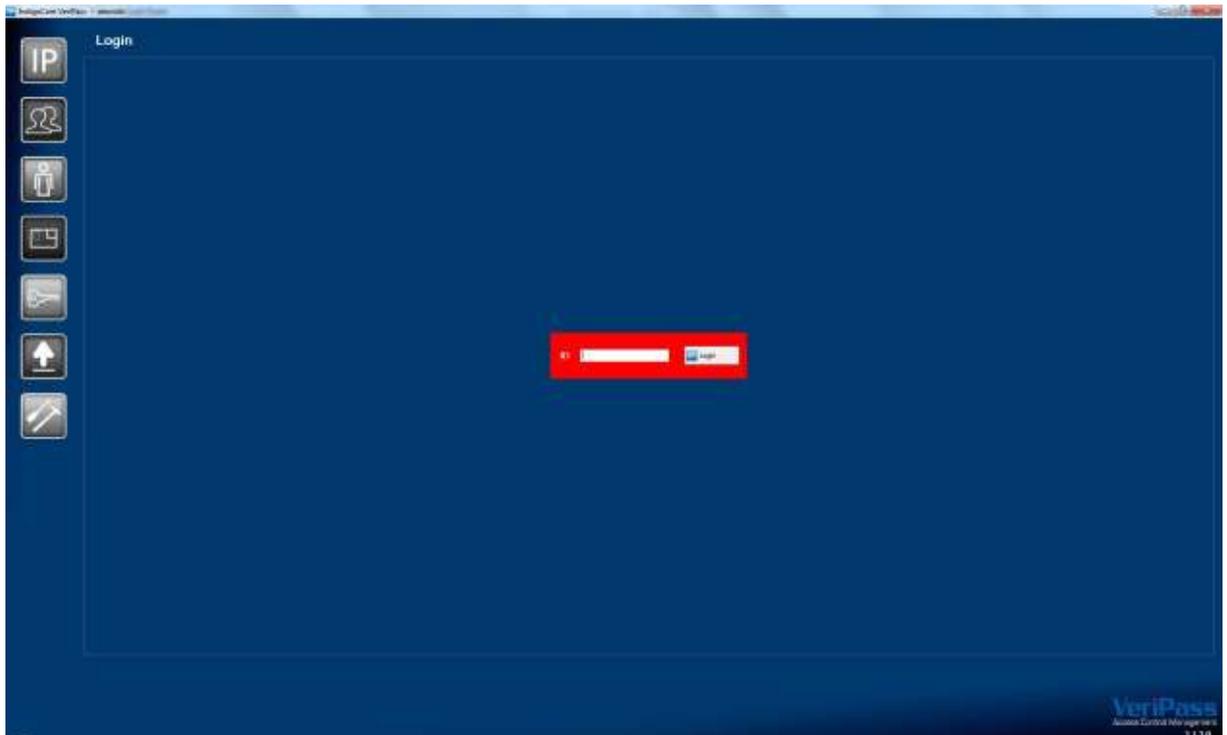

Before you can read or modify access rights, you have to login.

The default (installation) password is “icVeriPass”. For other user access, you can give administration rights to any user, so they can login with their appropriate iButton/RFID.


For each project it’s recommended to use the VeriPass on 1 dedicated PC, because not all settings can be retrieved from the access control units. Some settings are locally saved.

3. Menu

	Scan the network for Access Control units
	Manage groups
	Manage users
	Manage rooms
	Manage access rights
	Upload all settings to the system
	Settings screen

Click on the respective icon to open its settings.

3.1. Scan the network for access control units

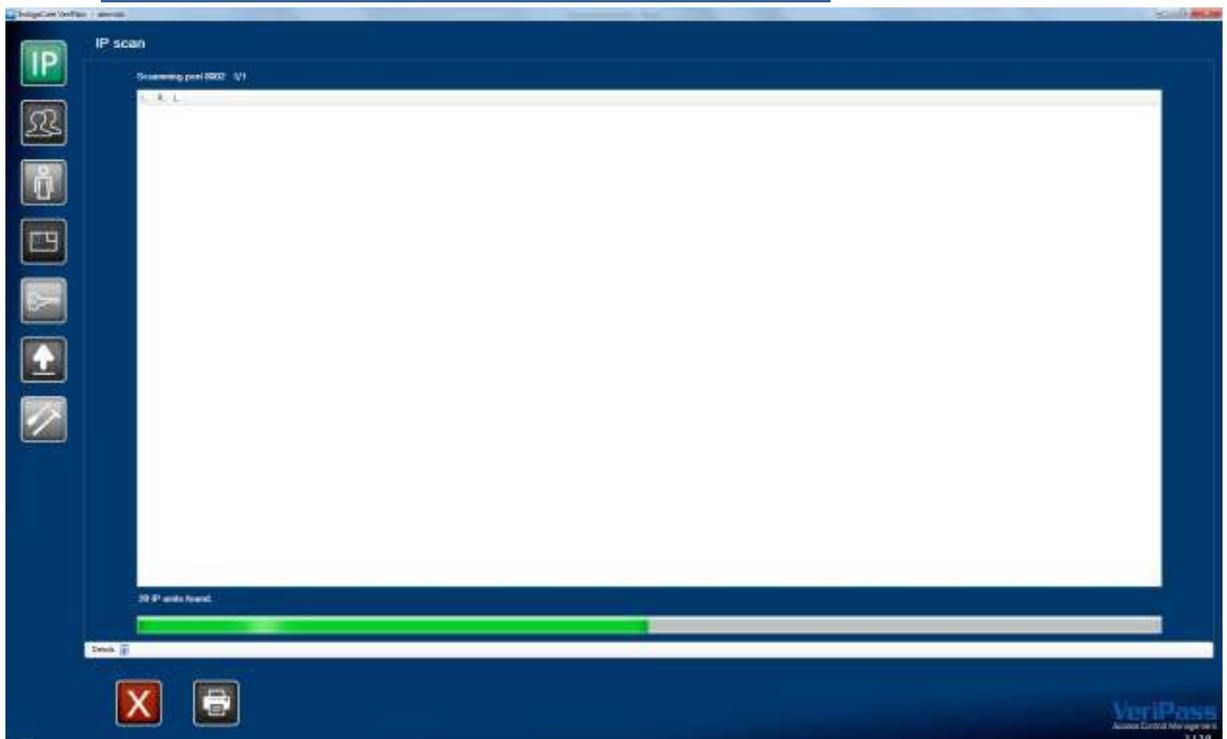

On startup, the application will automatically scan the network for access control units. It will show differences with the latest scan in a listview.




	No difference
	Not found
	New

Rescans can be made manually by pressing the rescan button.

	Cancel network scan
	Printable overview
	Scan again


3.2. Manage groups

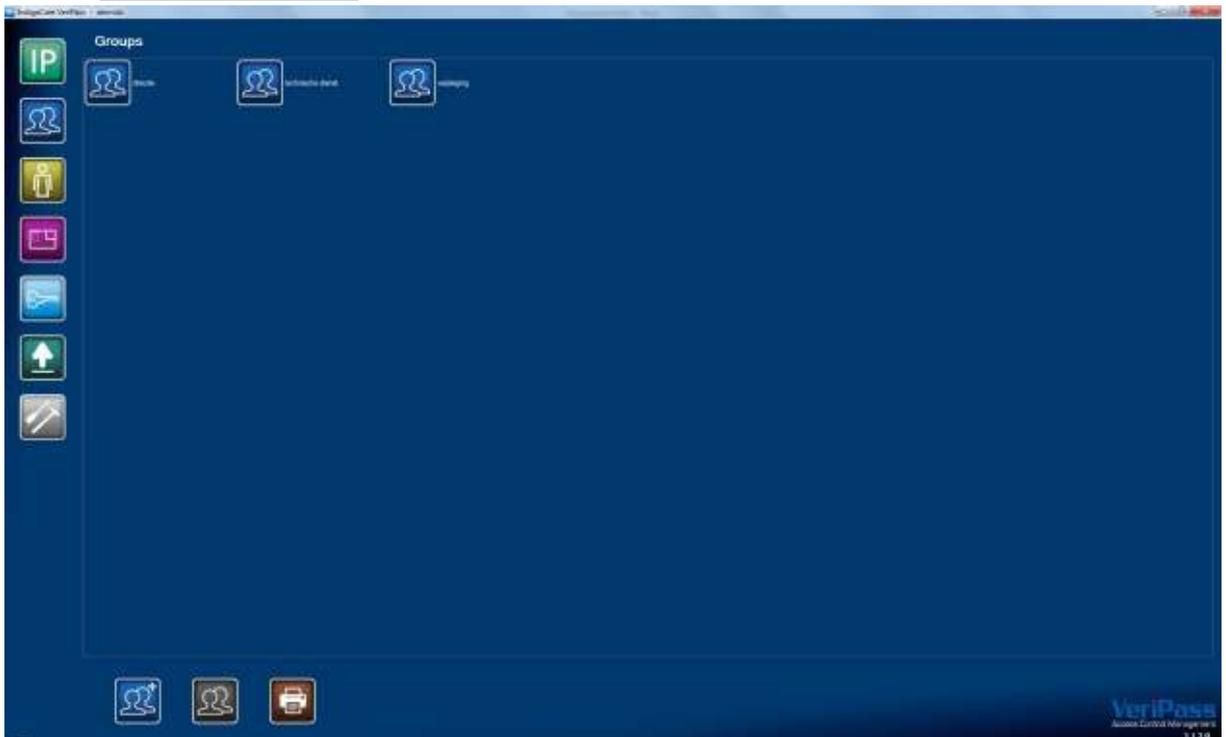

This section shows all existing groups. A group is a collection of users that have the same access rights. Groups can be changed, added or deleted. To edit a group, double click on its icon. To add a user to a group, go to the user menu, select the appropriate user and select the appropriate group. It is recommended to create your groups before you start to create your users.

	Add group
	Delete group
	Printable overview

Group settings

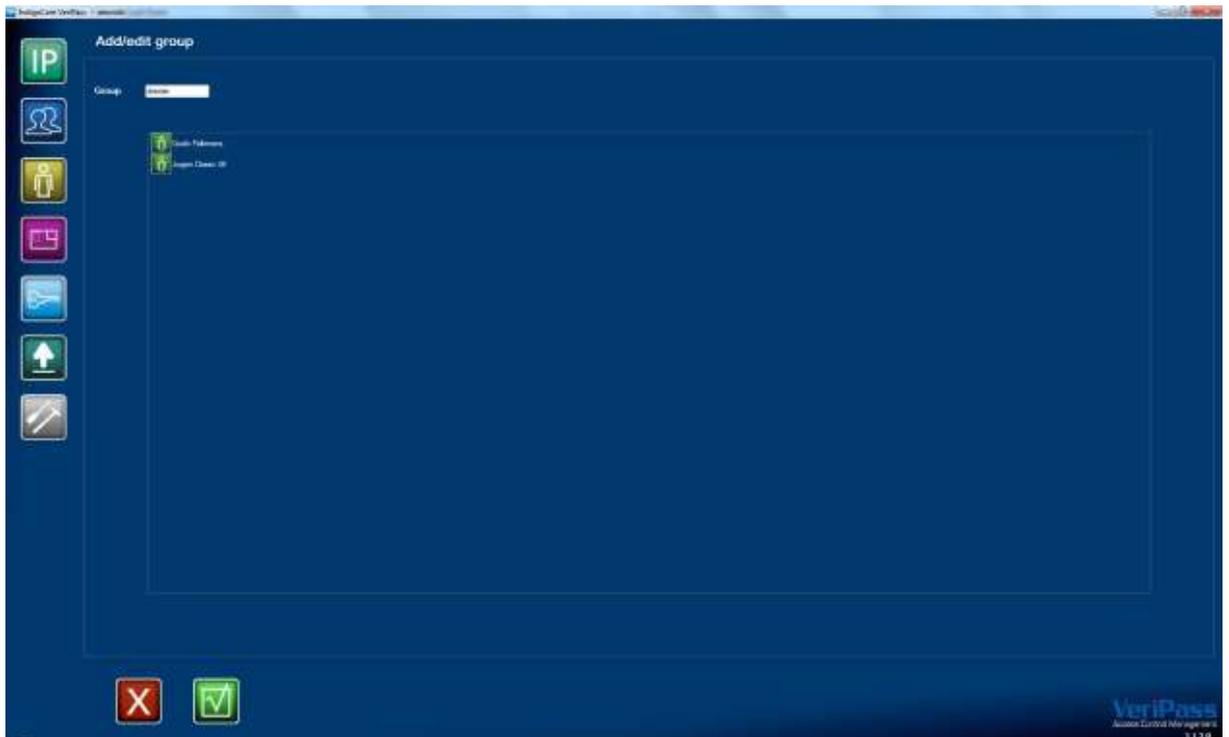

All group settings can be made here. An overview of its users is available.

	Cancel adjustments
	Save adjustments

3.3. Manage users



This section shows all existing users. Users can be changed, added or deleted. To edit a user, double click on its icon.

	Add user
	Delete user
	Printable overview
	Import user list from iLink
	Export user list to iLink
	Export to CSV

In the upper right corner is search field is added to find users or an iButton/RFID. Type a part of the name or scan your ID in this field to find the appropriate user.



User settings




All user settings can be made here. To add an ID, you can type your ID manually, or you can scan the ID with your iButton or RFID reader.

Select the Admin checkbox to grant admin rights to the user. A user with admin rights, can login in this software (with its badge) to change the settings.

Below the admin checkbox, 4 domotics rights can be granted for iCall Sip Touch units. The name for these domotics can be set in the setup section.

	Cancel adjustments
	Save adjustments


3.4. Manage rooms

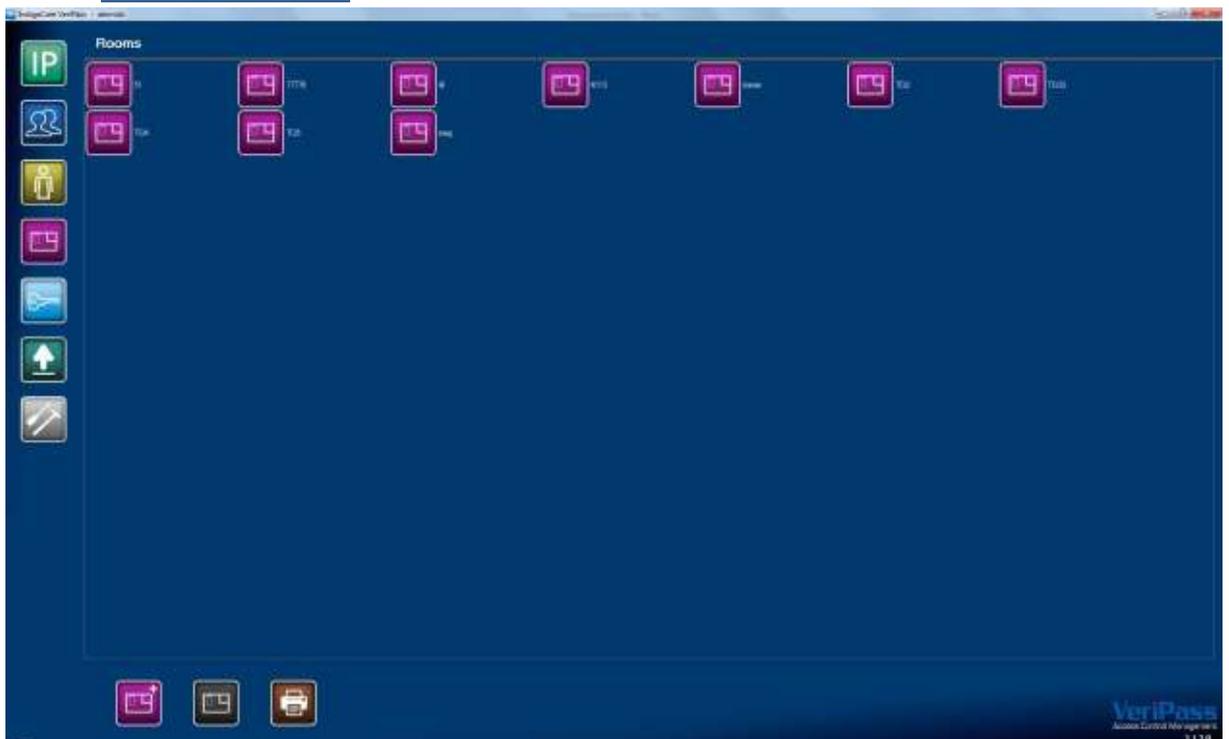

This section shows all existing rooms. Rooms can be changed, added or deleted. To edit a room, double click on its icon.

	Add room
	Delete room
	Printable overview


Room settings

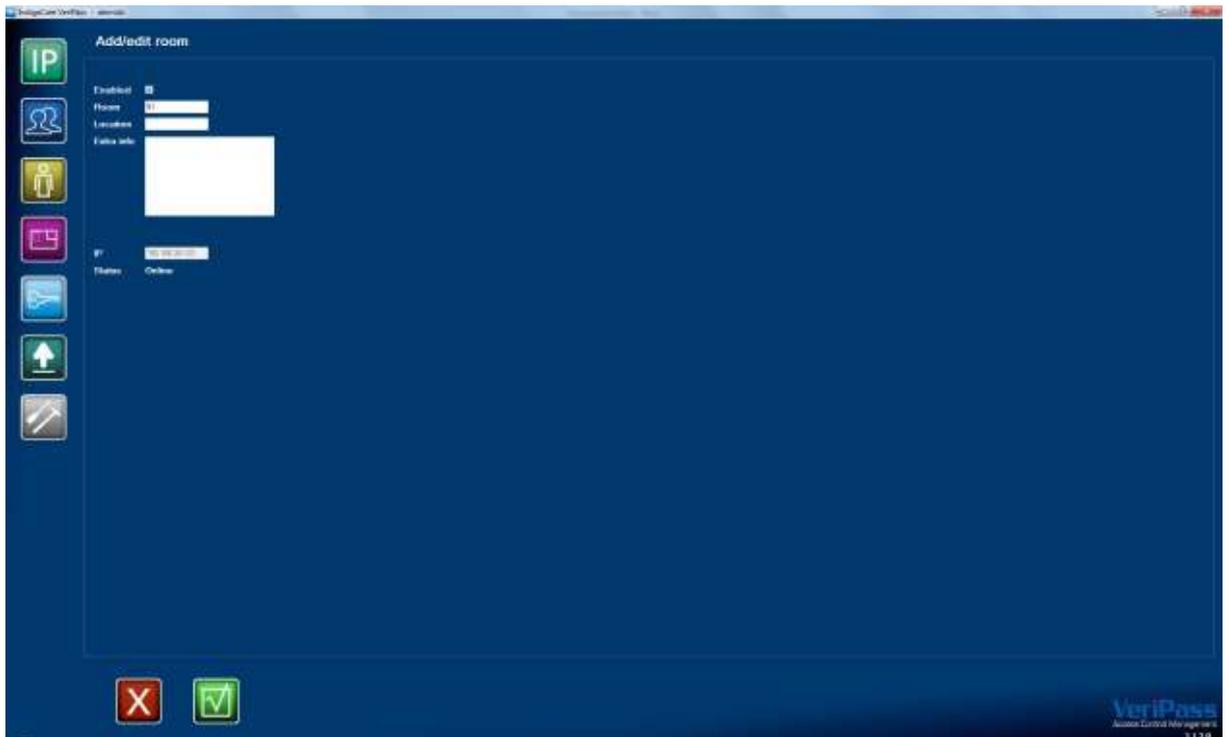

All room settings can be made here.

	Cancel adjustments
	Save adjustments


3.5. Manage access rights

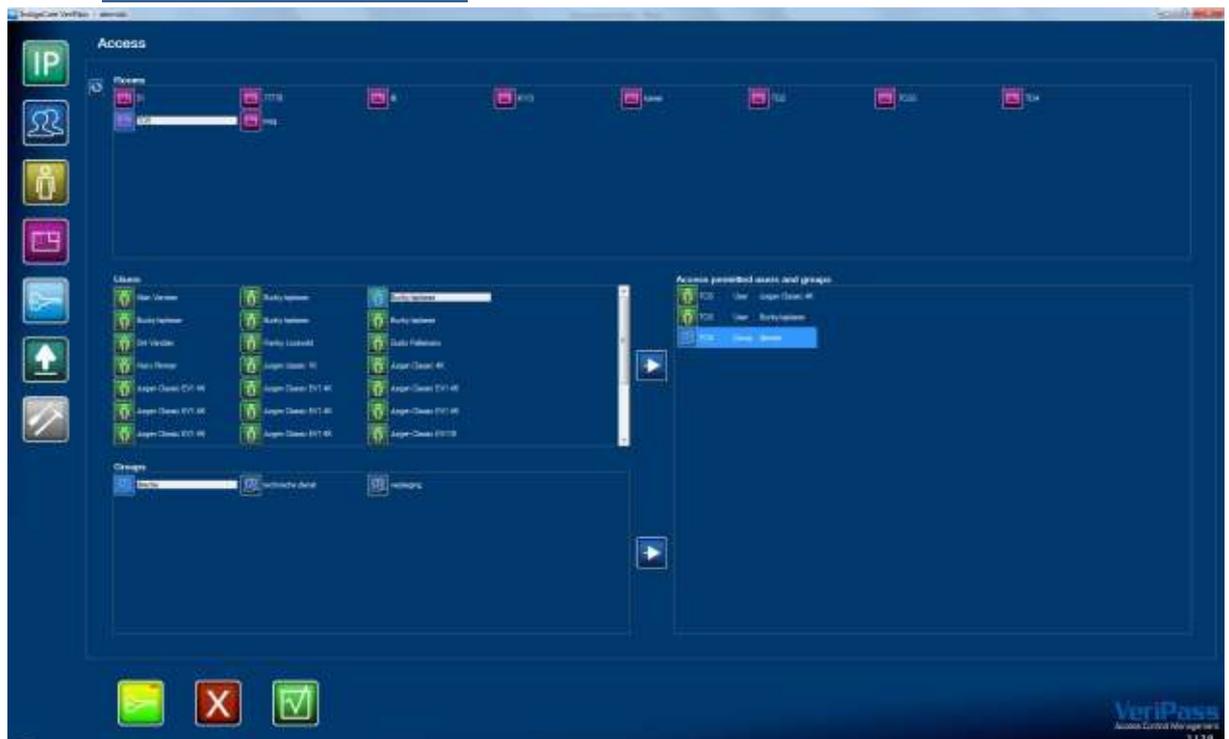

In this section all access rights can be granted or deleted. The new access rights will only be active if you upload the new settings to the IP units.

	Delete user/group access right
	Cancel adjustments
	Save adjustments
	Add access right


3.6. Upload all settings to the system

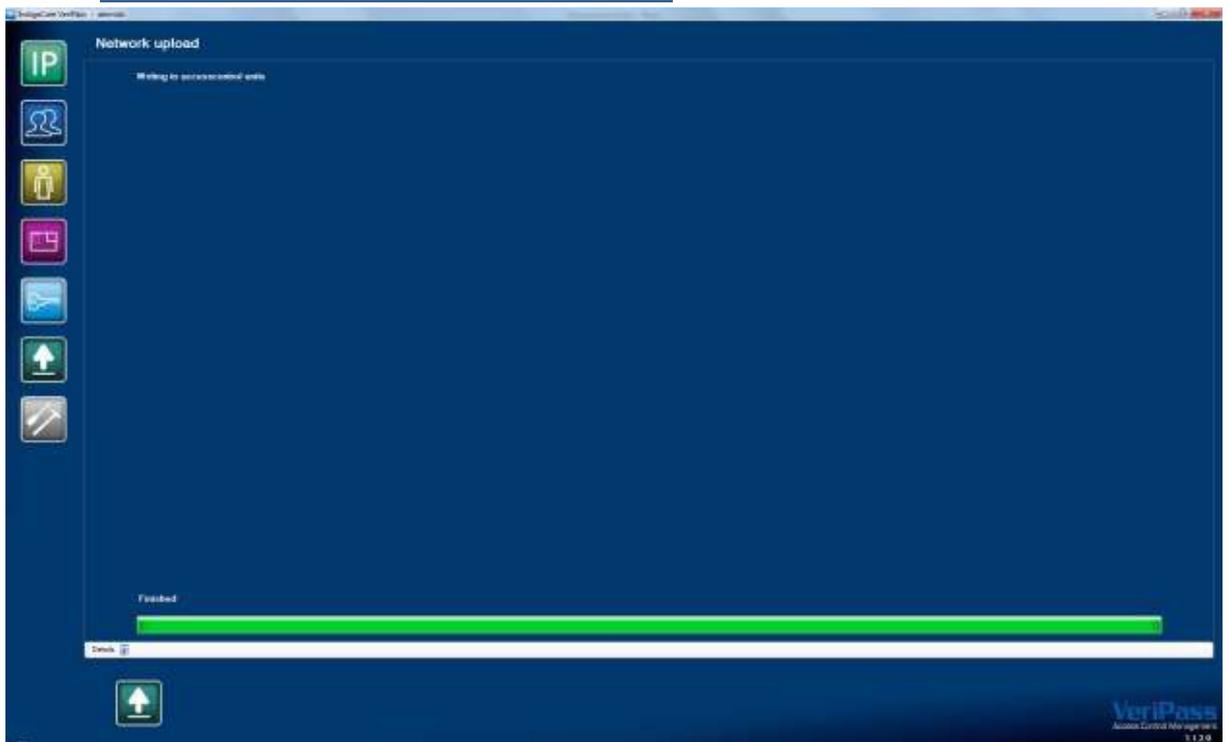

In this section all access rights can be sent to the IP units. The new access rights will be active when the upload is ready.

	Start upload
	Cancel upload

3.7. Settings screen

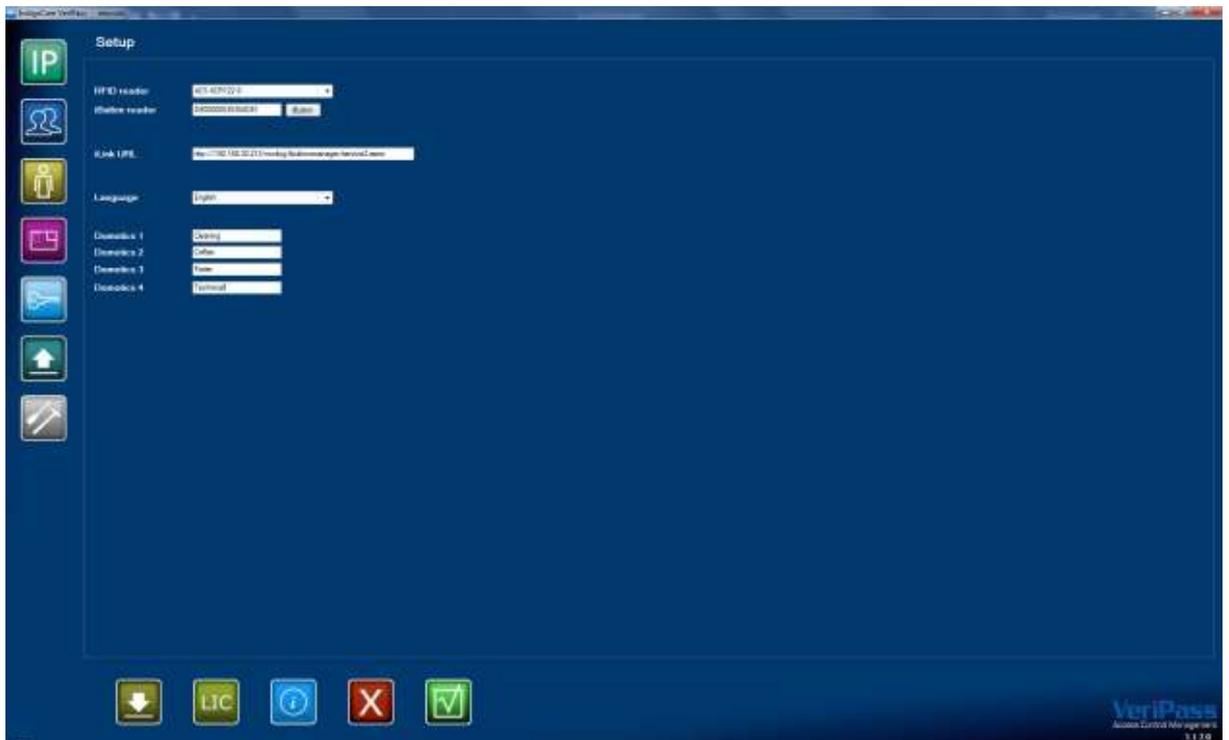

In this section basic settings can be saved.

	Backup
	License settings
	About the program
	Cancel adjustments
	Save adjustments